




Revo[®] 3-wheel

Revo[®] 4-wheel

Model number:	SC63	SC64
Weight capacity:	300 lbs.	300 lbs.
Maximum speed:	Up to 4.5 mph	Up to 4.5 mph
Ground clearance:	2.25"	2.25"
Turning radius:	38"	50"
Overall length:	43"	44.75"
Overall width:	21.5"	21.5"
Wheels		
Front:	2" x 8" solid	2" x 8" solid
Rear:	3" x 9" solid	3" x 9" solid
Anti-tips:	Rear	Rear
Suspension:	N/A	N/A
Transport weight		
Total w/o batteries:	106.5 lbs.	125 lbs.
Heaviest piece when disassembled:	48 lbs. (rear section)	51.5 lbs. (front section)
Drive train:	Sealed transaxle, 24 volt DC motor	
Dual braking system:	Electronic, regenerative and electromechanical (electronic, regenerative disc brakes)	
Power elevating seat:	N/A	
Std. Seating:	Compact foldable molded plastic	
Battery charger:	Off-board, 5 amp (standard)	
Per-charge range:	Up to 20 miles	
Battery requirements:	(2) 12V deep cycle (AGM or gel cell type recommended)	
Battery size:	U-1	
Warranty:	3-year limited warranty	

Colors:

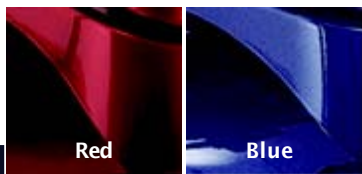

Red

Blue

Benefits

- Innovative interlocking drop-in battery packs for easy removal/installation
- Easily disassembles into 5 lightweight, easy-to-manage pieces
- 300 lbs. weight capacity
- Easy-to-use external freewheel lever
- State-of-the-art, microprocessor-based motor controller
- Seat offers height adjustment, adjustable 360° seat rotation, and flip-up armrests for optimal user comfort
- Adjustable tiller for maximum comfort
- Headlight and taillight for added safety
- Easily accessible charger port and circuit breaker

Options

- Forearm crutch holder
- Lap belt
- Cane/crutch holder
- Wishbone crutch holder
- Oxygen tank holder
- Cup holder
- Rear basket
- Safety flag
- Walker holder
- Weather cover
- Saddlebag
- Quad cane holder
- Cane, clip attachment



- Backpacker
- Double Entry
- Full Platform



Backpacker



Double entry

These lifts are compatible with the Revo[®].

The information contained herein is correct at the time of publication; we reserve the right to alter specifications without prior notice. Speed & range vary with user weight, terrain type, battery charge, battery condition and tire pressure.

800-800-8586 (US) • 888-570-1113 (Canada) • www.pridemobility.com



REVO®

THE REVO PROVIDES SLEEK STYLING
AND SUPER PORTABILITY.

REVO®

The Revo® delivers “full-sized” scooter style, performance, and features in a compact, lightweight design that easily disassembles for transport while remaining highly affordable.



Pride[®]
Mobility Products Corp.

800-800-8586 (US) · 888-570-1113 (Canada) · www.pridemobility.com

REVOSS01-07